

ME 100



The ideal extractor for light industry and laboratory environments.

With its optimal design, \emptyset 4 inch Movex ME has a very low pressure drop, which provides many valuable benefits.

- Low pressure drop saves energy.
- Air flow noise is reduced.
- Lower pressure drop is achieved without selecting a larger diameter extractor.
- Lower pressure drop allows the ME to be combined with additional extraction systems.

To further facilitate maneuvering of the extractor, the models 1650 and 1900 are equipped with a pulling gas spring as standard and the models 2100 and 2650 with two pulling gas springs.

An easy-to-grip handle facilitates the maneuverability of the extractor.

Unique design and stable mounting brackets make the Movex ME your best choice Support for designing an effective system can be found on page 5, and at www.movexinc.com where you will find our design tool and CAD drawings.

The Movex range also includes fans, accessories, automatic controls and filters suitable for local extraction.

LOCAL EXTRACTION Pure advantages

Always choose a low pressure drop

Lowest possible pressure drop is a quality aspect that always should be considered.


With its uniquely designed joint construction, Movex ME combines maximum flexibility with low pressure drop. The air passes through the joints without creating unnecessary turbulence, thus producing an energy-saving low pressure drop and a quieter working environment.



Recommended air flow

The recommended air flow for a Ø4 arm is 120-180 cfm, See table and diagram.

Activity	Air flow
Laboratories	120-180 cfm
Light industry	180 cfm

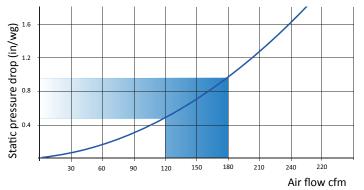

Static pressure drop is measured in accordance with ISO standard 5167-1.

Optimal capture

For optimum benefit from the local extractor, it is important to use the flexibility of the extractor to get as close to the contaminant as possible. A good rule of thumb would be a distance of 2–3 times the diameter of the extractor tube. At the recommended air flow, the extractor will provide high efficiency even if disturbances are generated in the surroundings.

Unique benefits

The Movex ME's joints have a patented friction design that, combined with the large joint diameter and single grip handle, provide a secure, position-stable arm with smooth adjustments. All without the need to apply excessive force or use tools on the adjusting knob.

Joints with reinforced ends and ball bearings moderate the friction and allow the arm to be moved up and down while maintaining stability and function.

Handle for easy adjustment

A steady and easily accessible handle, that provides easy adjustment, is fitted as a standard on all models of Movex ME \emptyset 4 inch arm.



One arm. All options.

Movex ME has a complete range of accessories to suit every situation, enabling you to create the optimal extractor for the evacuation of hazardous airborne gases and particulates.



Standard version

Suitable for evacuating most types of airborne contaminants, e.g. in laboratories, schools, hospitals, the pharmaceutical industry, nail salons and light industrial applications.



PP version

Used primarily for evacuating very corrosive contaminants in high concentrations, e.g. in certain laboratories and in the pharmaceutical and chemical industries.




ATEX version Suitable for evacuating airborne contaminants where there is a requirement for an ATEX-classified environment, e.g. in laboratories, the chemical and petrochemical industries, gas distribution, and the paint and pharmaceutical industries.



MET for ceiling and wall mounting, 3 joints

Standard	Size (inches)			Gas spring	Weight
	Α	В	øс	pcs.	(lb)
MET 1150-100	18	14	4	0	10.8
MET 1350-100	22	18	4	0	11.9
MET 1650-100	30	22	4	1	13
MET 1900-100	40	22	4	1	14.1
MET 2100-100	40	30	4	2	15.2
MET 2650-100	52	40	4	2	16.3
РР	Si	ze(inch	es)	Gas spring	Weight
РР	Si: A	ze(inch B	es) øc	Gas spring pcs.	Weight (Ib)
PP MET 1150-100PP					•
	Α	В	øc	pcs.	(lb)
MET 1150-100PP	А 18	в 14	øс 4	pcs.	(lb) 10.8
MET 1150-100PP MET 1350-100PP	A 18 22	в 14 18	øс 4 4	рся. О О	(Ib) 10.8 11.9
MET 1150-100PP MET 1350-100PP MET 1650-100PP	A 18 22 30	в 14 18 22	øc 4 4 4	0 0 1	(lb) 10.8 11.9 13

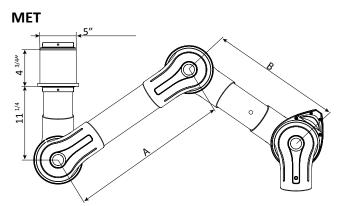


Add MTI bracket for ceiling mounting. Add MVK bracket for wall mounting.

ATEX	Size (inches)			Gas spring	Weight
	Α	В	øс	pcs.	(lb)
MET 1150-100EX	18	14	4	0	10.8
MET 1350-100EX	22	18	4	0	11.9
MET 1650-100EX	30	22	4	1	13
MET 1900-100EX	40	22	4	1	14.1
MET 2100-100EX	40	30	4	2	15.2
MET 2650-100EX	52	40	4	2	16.3

	Total length of arm
	(inches)
ME 1150-100	46"
ME 1350-100	54"
ME 1650-100	66"
ME 1900-100	76″
ME 2100-100	84"
ME 2650-100	106″

Pressure drop

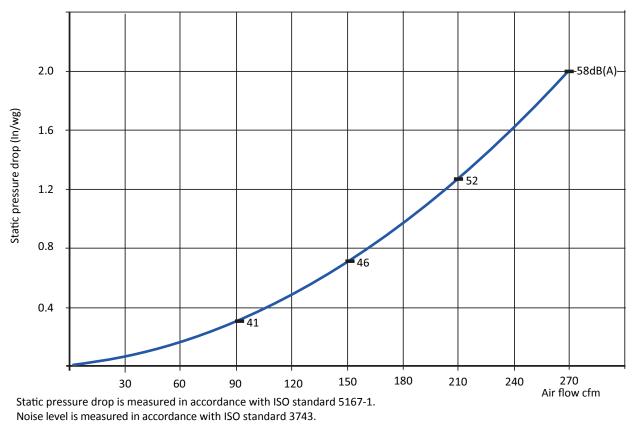

Indicated sound level refers to sound pressure level.

MEB for table mounting, 3 joints

Standard	Siz	e (inche	es)	Gas spring	Weight	N
	Α	В	øс	pcs.	(lb)	
MEB 1150-100	18	14	4	0	10.8	
MEB 1350-100	22	18	4	0	11.9	
MEB 1650-100	30	22	4	1	13	
MEB 1900-100	40	22	4	1	14.1	
РР	Siz	e (inche	es)	Gas spring	Weight	
	А	В	øс	pcs.	(lb)	
MEB 1150-100PP	18	14	4	0	10.8	
MEB 1350-100PP	22	18	4	0	11.9	1/1
MEB 1650-100PP	30	22	4	1	13	,
MEB 1900-100PP	40	22	4	1	14.1	
ATEX	Siz	e (inche	es)	Gas spring	Weight	
	А	В	øc	pcs.	(lb)	

14

18

22

22

4

4

4

4

0

0

1


1

10.8

11.9

13

14.1

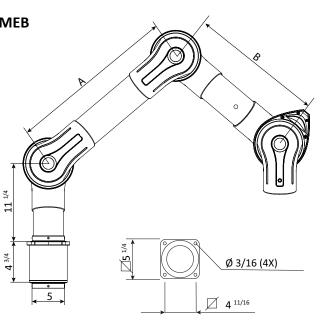

Reach at recommended installation height

18

22

30

40

The following installation heights and sideways displacement relative to the work area are recommended for optimal extraction:

Recommended i	installation	height
---------------	--------------	--------

MEB 1150-100EX

MEB 1350-100EX

MEB 1650-100EX


MEB 1900-100EX

Designation	H (inches)
MET 1150-100	68-76
MET 1350-100	76-84
MET 1650-100	80-88
MET 1900-100	88-96
MET 2100-100	92-100
MET 2650-100	92-100

Recommended side displacement


radius, relative to work area

Designation	S (inches)
MET 1150-100	12-24
MET 1350-100	16-28
MET 1650-100	20-32
MET 1900-100	28-32
MET 2100-100	28-36
MET 2650-100	52-36

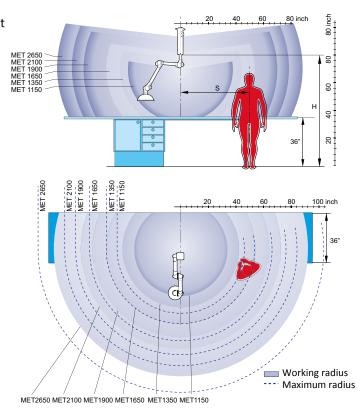

Hoods



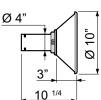

METAL HOOD

The metal hood is used when working in corrosive environments and for capturing hot gasses and dust splatter. Metal hoods can be fitted with work lighting. Temp. range: -5°F to +176°F

Standard	Variants	Weight (oz)
MEM 251-100	PP, EX	18

Material

Standard/PP: Powder-coated aluminum ЕΧ Powder-coated aluminum





4

Ø


The clear dome hood is suitable for lighter gasses with a wider dispersal of contaminants without blocking the user's vision.

-5°F to +176°F Temp. range:

Standard	Variants	Weight (oz)
MEK 351-100	PP,EX	21.5

Material	
Standard:	PMMA
PP	Polypropylene, transparent
EX	PEEL black



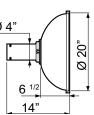


DOME HOOD

The larger clear dome hood is also suitable for lighter gasses with an even wider dispersal of contaminants, still without blocking the user's vision. Temp. range: -5°F to +176°F

Standard	Variants	Weight (oz)
MEK 500-100	PP,EX	26
otanidarai	PMMA Polypropylene, c	ppaque

Standard:	PMMA
PP	Polypropylene, opaque
EX	PEEL black

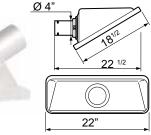


SQUARE HOOD

The square hood is suitable for placing above gases with a high lift, or adjacent to the work surface for contaminants with no lift or low lift - all this without interfering with the work. Temp. range: -5°F to +176°F

Standard	Variants	Weight (oz)
MESH 500-100		40

Material Standard: PETG



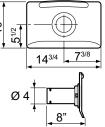

FLAT SCREEN HOOD

The flat screen hood is designed to maximise the working area without obscuring the object from the user. The flat screen hood gives the best suction effect for table and bench tasks. Temp. range: -5°F to +176°F

	(oz)
PP,ES, EX	22
PETG	
olypropylene	
,	ETG



PROTECTIVE GRILL

The protective grill (mounted in the first joint) keeps objects out of the extractor. Temp. range: -5°F to +176°F


Standard	Variants	Weight (oz)
MSG-100	EX	0.40

Brackets




All Movex laboratory extractors have as standard a full swivel that allows 360° of rotation without the need to add special sleeve couplings.

Both ceiling and wall brackets have a special squareshaped profile in anodized aluminum to provide a stylish and stable installation. This aluminum profile also allows both the wall and ceiling brackets to be custom tailored at the job site.

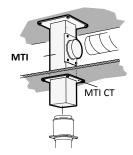

The MTI ceiling bracket

The ceiling bracket functions as a simple and stable duct for outgoing air, avoiding the need for expensive ducting and additional holes through false ceilings. **On request, the MTI can be supplied in lengths exceeding 80 inches.**

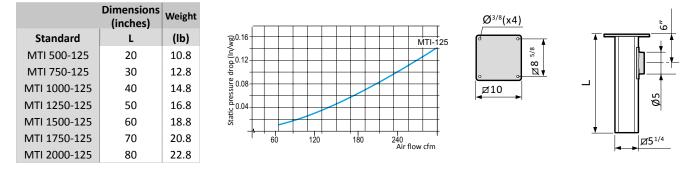


The MTF ceiling bracket

Ceiling bracket, for fitting through beams. The attachment plate is adjustable for the entire length of the aluminum profile. If required, the aluminum profile can be cut during fitting.

	Dimensions (inches)	Weight
Standard	L	(lb)
MTF-125	40	9.3

As well as the standard design, the MTI/MTF is available in an ATEX (EX) version. The ceiling brackets can be supplied with an epoxy-coated exterior in all lengths. For aggressive environments, we recommend epoxy coating on the interior and exterior.

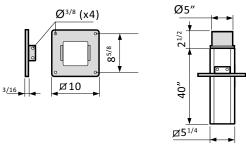


The MTI CT escutcheon plate

Escutcheon plate, used with the MTI ceiling brackets for stabilization and to cover the rough cut in the false ceiling.

	Weight
Standard	(oz)
MTI CT-125	4.4

As well as the standard design, the escutcheon plate is available in an ATEX (EX) version.

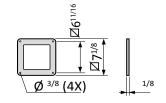



MVK wall bracket

Wall brackets can be special ordered in custom horizontal and vertical lengths.

	Weight
Standard	(lb)
MVK-125	3.4

As well as the standard design, the wall bracket is available in an ATEX (EX) version.

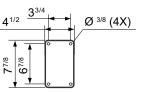

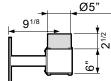



Material description

Friction joints

Ball bearing-equipped adjustable friction joints in polypropylene (PP), with guide ring in low friction-treated rubber.Support springs and other component parts in zinc-plated steel or stainless steel.

Tubes

Made from thin-walled anodized aluminum or from polypropylene. Air-tight damper supplied as standard.

ME Standard

The standard ME version has polypropylene joints and anodized aluminum tubes.

The standard ME version is suitable for evacuating most types of airborne contaminants, e.g. in laboratories, schools, hospitals, the pharmaceutical industry, nail salons and light industrial applications.

ME PP

Polypropylene joints and tubes version. All metallic parts that are in contact with the air flow are made of stainless steel.

The PP version of the ME is used primarily for evacuating very corrosive contaminants in high concentrations, e.g. in certain laboratories and in the pharmaceutical and chemical industries. When using a PP extractor fitted to a ceiling, we recommend that you order the MTI ceiling bracket with an internal epoxy coating.



Conductive polypropylene joints and tubes. All metallic parts that are in contact with the air flow are made of stainless steel. Static electricity is diverted to a separate grounded connection. All steel supporting parts are lined in a conductive powder coating. The product meets the requirements of category 2 of the ATEX directive (94/9/EC) for gases and dust.

The ATEX version of the ME is suitable for evacuating airborne contaminants where there is a requirement for ATEX-classified products, e.g. laboratories, the chemical and petrochemical industries, gas distribution, and the paint and pharmaceutical industries.

Delivery

Ceiling/Supplied assembled, excluding hood. TheWallceiling or wall brackets should be ordered separately.MET

Table-Supplied assembled, with attachment plate for tableMEBfitting, excluding hood.The MBF flexible table bracket should
be ordered separately.



